

Mercodia Porcine Insulin ELISA





Directions for Use

10-1200-01

REAGENTS FOR 96 DETERMINATIONS

Manufactured by
Mercodia AB,
Sylveniusgatan 8A,
SE-754 50 Uppsala,
Sweden

EXPLANATION OF SYMBOLS USED ON LABELS

 $\Sigma = 96$	Reagents for 96 determinations
	Expiry date
	Store between 2-8°C
	Lot No.

INTENDED USE

Mercodia Porcine Insulin ELISA provides a method for the quantitative determination of insulin in porcine serum and plasma.

SUMMARY AND EXPLANATION OF THE TEST

Insulin is the principal hormone responsible for the control of glucose metabolism. It is synthesised in the β -cells of the islets of Langerhans as the precursor, proinsulin, which is processed to form C-peptide and insulin. Both are secreted in equimolar amounts into the portal circulation. The mature insulin molecule comprises two polypeptide chains, the A chain and the B chain. The two chains are linked together by two inter-chain disulphide bridges. There is also an intra-chain disulphide bridge in the A chain.

Secretion of insulin is mainly controlled by plasma glucose concentration, and the hormone has a number of important metabolic actions. Its principal function is to control the uptake and utilisation of glucose in peripheral tissues via the glucose transporter. This and other hypoglycaemic activities, such as the inhibition of hepatic gluconeogenesis and glycogenolysis are counteracted by the hyperglycaemic hormones including glucagon, epinephrine (adrenaline), growth hormone and cortisol.

Islet cell transplantation has long been considered as a potential cure for type 1 diabetes. The shortage of human donors and difficulty in isolating purified islets from adult human pancreas has drawn the attention to porcine pancreatic islets (1-4). Pig and human insulin are structurally similar, and the regulation of insulin secretion in pigs resembles humans. Parameters such as insulin secretion are used to investigate the viability of the transplanted islets (4-6).

PRINCIPLE OF THE PROCEDURE

Mercodia Porcine Insulin ELISA is a solid phase two-site enzyme immunoassay. It is based on the direct sandwich technique in which two monoclonal antibodies are directed against separate antigenic determinants on the insulin molecule. During incubation, insulin in the sample reacts with peroxidase-conjugated anti-insulin antibodies and anti-insulin antibodies bound to the microplate. After a simple washing step that removes unbound enzyme labelled antibody, the bound conjugate is detected by reaction with 3,3',5,5'-tetramethylbenzidine (TMB). The reaction is stopped by the addition of acid, giving a colorimetric endpoint that can be read spectrophotometrically.

WARNINGS AND PRECAUTIONS

- The content of this kit and their residues must not be allowed to come into contact with ruminating animal or swine.
- The Stop Solution in this kit contains 0.5 M H_2SO_4 . Follow routine precautions for handling hazardous chemicals.
- All samples should be handled as capable of transmitting infections.
- Not for internal or external use in humans or animals.
- Each well can only be used once.

MATERIAL REQUIRED BUT NOT PROVIDED

- Pipettes with appropriate volumes (repeating pipettes preferred for addition of enzyme conjugate 1X solution, Substrate TMB and Stop Solution)
- Tubes, beakers and cylinders for reagent preparation
- Redistilled water
- Magnetic stirrer
- Vortex mixer
- Microplate reader with 450 nm filter
- Microplate shaker (700–900 cycles per minute, orbital movement)
- Microplate washing device with overflow function (recommended but not required)

REAGENTS

Each Mercodia Porcine Insulin ELISA (10-1200-01) kit contains reagents for 96 wells, sufficient for 42 samples and one calibrator curve in duplicate. For larger series of assays, use pooled reagents from packages bearing identical lot numbers. The expiry date for the complete kit is stated on the outer label. The recommended storage temperature is 2–8°C.

Coated Plate Mouse monoclonal anti-insulin For unused microplate strips, reseal the bag using adhesive tape, store at 2–8°C and use within 8 weeks.	1 plate	96 wells 8-well strips	Ready for Use
Calibrators 1, 2, 3, 4, 5 Porcine insulin Color coded yellow Concentration stated on vial label	5 vials	1000 µL	Ready for Use
Calibrator 0 Color coded yellow	1 vial	5 mL	Ready for Use
Enzyme Conjugate 11X Peroxidase conjugated mouse monoclonal anti-insulin	1 vial	1.3 mL	Preparation, see below
Enzyme Conjugate Buffer Color coded blue	1 vial	13 mL	Ready for use
Wash Buffer 21X Storage after dilution: 2–8°C for 8 weeks	1 bottle	50 mL	Dilute with 1000 mL redistilled water to make wash buffer 1X solution.
Substrate TMB Colorless solution <i>Note! Light sensitive!</i>	1 bottle	22 mL	Ready for Use
Stop Solution 0.5 M H ₂ SO ₄	1 vial	7 mL	Ready for Use

Preparation of enzyme conjugate 1X solution

Prepare the needed volume of enzyme conjugate 1X solution by dilution of Enzyme Conjugate 11X, (1+10) in Enzyme Conjugate Buffer according to the table below. Mix gently.

When preparing enzyme conjugate 1X solution for the whole plate or if the reagents are to be used within 2 weeks, pour all of the Enzyme Conjugate Buffer into the Enzyme Conjugate 11X vial.

Number of strips	Enzyme Conjugate 11X	Enzyme Conjugate Buffer
12 strips	1 vial	1 vial
8 strips	700 µL	7 mL
4 strips	400 µL	4 mL

Storage after dilution: 2–8°C for 2 weeks.

SPECIMEN COLLECTION AND HANDLING

Serum

Collect blood by venipuncture, allow to clot and separate the serum by centrifugation. Samples can be stored at 2–8°C up to 24 hours. For longer periods, store samples at –20°C. Avoid repeated freezing and thawing.

Plasma

Collect blood by venipuncture into tubes containing heparin, citrate or EDTA as anticoagulant, and separate the plasma fraction. Samples can be stored at 2–8°C up to 24 hours. For longer periods store samples at –20°C. Avoid repeated freezing and thawing.

PREPARATION OF SAMPLES

No dilution is normally required for serum or plasma. All samples containing porcine insulin above the highest calibrator should be diluted with Calibrator 0 or with Mercodia Diabetes Sample Buffer (10-1195-01).

TEST PROCEDURE

All reagents and samples must be brought to room temperature before use. Perform each determination in duplicate for calibrators and samples. Prepare a calibrator curve for each assay run.

1. Prepare enzyme conjugate 1X solution (according to the table on previous page) and wash buffer 1X solution.
2. Prepare sufficient microplate wells to accommodate Calibrators, controls and samples in duplicate.
3. Pipette 25 μ L each of Calibrators, controls and samples into appropriate wells.
4. Add 100 μ L of enzyme conjugate 1X solution into each well.
5. Incubate on a plate shaker (700-900 rpm) for 2 hours at room temperature (18-25°C).
6. Wash 6 times with 700 μ L wash buffer 1X solution per well using an automatic plate washer with overflow-wash function, after final wash, invert and tap the plate firmly against absorbent paper. Do not include soak step in washing procedure.
Or manually,
Discard the reaction volume by inverting the microplate over a sink. Add 350 μ L wash solution to each well. Discard the wash solution, tap firmly several times against absorbent paper to remove excess liquid. Repeat 5 times. Avoid prolonged soaking during washing procedure.
7. Add 200 μ L Substrate TMB into each well.
8. Incubate on the bench for 15 minutes at room temperature (18-25°C).
9. Add 50 μ L Stop Solution to each well.
Place the plate on the shaker for approximately 5 seconds to ensure mixing.
10. Read optical density at 450 nm and calculate results.
Read within 30 minutes.

Note! To prevent contamination between the conjugate and substrate, separate pipettes are recommended.

INTERNAL QUALITY CONTROL

Commercial controls such as Mercodia Animal Insulin Control (10-1221-01) and/or internal serum pools with low, intermediate and high porcine insulin concentrations should routinely be assayed as samples, and results charted from day to day. It is good laboratory practice to record the following data for each assay: kit lot number, dilution and/or reconstitution dates of kit components, OD values for the blank, Calibrators and controls.

Laboratories should follow government regulations or accreditation requirements for quality control frequency.

CALCULATION OF RESULTS

Computerized calculation

The concentration of porcine insulin is obtained by computerized data reduction of the absorbance for the Calibrators, except Calibrator 0, versus the concentration using cubic spline regression.

Manual calculation

1. Plot the absorbance values obtained for the Calibrators, except Calibrator 0, against the porcine insulin concentration on a log-log paper and construct a calibrator curve.
2. Read the concentration of the samples from the calibrator curve.


Example of results

Wells	Identity	A ₄₅₀	Mean conc. µg/L
1A-B	Calibrator 0	0.054/0.053	
1C-D	Calibrator 1*	0.090/0.086	
1E-F	Calibrator 2*	0.135/0.133	
1G-H	Calibrator 3*	0.319/0.329	
2A-B	Calibrator 4*	0.907/0.880	
2C-D	Calibrator 5*	2.405/2.304	
2E-F	Sample 1	0.184/0.185	0.079
2G-H	Sample 2	0.571/0.560	0.290
3A-B	Sample 3	1.597/1.543	0.977

*Exact concentration stated on vial label.

Example of calibrator curve

A typical calibrator curve is shown here. Do not use this curve to determine actual assay results.

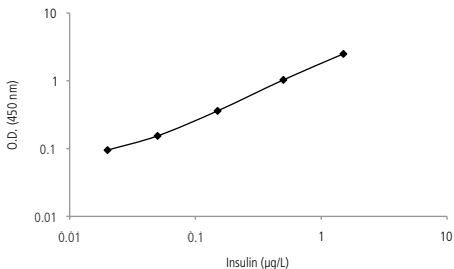

LIMITATIONS OF THE PROCEDURE

Grossly lipemic, icteric or hemolyzed samples do not interfere in the assay. However, hemolysis in serum and plasma samples may result in a degradation of insulin which could give falsely low values and contributes to higher inter assay variation. The degradation is dependent on time, temperature and the hemoglobin concentration. Keep hemolyzed samples cold or on ice to prevent the insulin degradation.

EXPECTED VALUES

Good practice dictates that each laboratory establishes its own expected range of values.

PERFORMANCE CHARACTERISTICS

Detection limit

Detection limit is defined as the Capability of Detection according to ISO11843-Part 1. Capability of Detection should be seen as part of a method validation, rather than the lowest concentration that can be measured.

The detection limit is 0.01 (µg/L) determined with the methodology described in ISO11843-Part 4.

Concentration of samples with absorbance below Calibrator 1 should not be calculated, instead expressed as less or equal to \leq the concentration indicated on the vial for Calibrator 1.

Recovery

Recovery upon addition is 100-111 % (mean 105 %).

Recovery upon dilution is 90-112 % (mean 97 %).

Hook effect

Samples with a concentration of up to 1000 µg/L can be measured without giving falsely low results.

Precision

Each sample was analyzed in 4 replicates on 13 different occasions.

Sample	Mean value (µg/L)	Coefficient of variation		
		within assay %	between assay %	total assay %
1	0.062	3.5	3.9	4.3
2	0.127	3.8	2.6	3.2
3	0.239	3.1	1.5	2.2

Specificity

The following cross reaction have been found:

Porcine C-peptide	< 0.001 %
Porcine Proinsulin	< 0.2 %
Human C-peptide	< 0.01 %
Human Insulin	22 %
Human Proinsulin	< 0.1 %
Humalog	< 0.0000006 %
Lantus	5.4 %
Levemir	< 0.0000002 %
NovoRapid	0.7 %
Equine Insulin	100 %
Canine Insulin	100 %

CALIBRATION

Mercodia Porcine Insulin ELISA is calibrated against an inhouse reference preparation of porcine insulin.

WARRANTY

The performance data presented here was obtained using the procedure indicated. Any change or modification in the procedure not recommended by Mercodia AB may affect the results, in which event Mercodia AB disclaims all warranties expressed, implied or statutory, including the implied warranty of merchantability and fitness for use.

Mercodia AB and its authorized distributors, in such event, shall not be liable for damages indirect or consequential.

REFERENCES

1. Bogdani M, Suenens K, Bock T, Pipeleers-Marichal M, In't Veld P, Pipeleers D (2005) Growth and Functional Maturation of beta-cells in Implants of Endocrine Cells Purified From Prenatal Porcine Pancreas. *Diabetes* 54:3387-3394
2. Dufrane D, D'hoore W, Goebbles RM, Saliez A, Guiot Y, Gianello P (2006) Parameters favouring successful adult pig islet isolations for xenotransplantation in pig-to-primate models. *Xenotransplantation* 13:204-214
3. Groth CG, Tibell A, Wennberg L, Bennet W, Lundgren T, Rydgård KJ, Lundin S, Lindeborg E, Korsgren O (2000) Clinical aspects and perspectives in islet xenotransplantation. *J Hepatobiliary Pancreat Surg* 7:364-369
4. Korsgren O, Jansson L, Eizirik D, Andersson A (1991) Functional and morphological differentiation of fetal porcine islet-like cell clusters after transplantation into nude mice. *Diabetologia* 34:379-386
5. Martinez-Puig D, Mourot J, Ferchaud-Roucher V, Anguita M, Garcia F, Krempf M, Pérez, JF (2006) Consumption of resistant starch decreases lipogenesis in adipose tissues but not in muscular tissues of growing pigs. *Livestock Science* 99:237-247
6. Ramsay TG, Bush JA, McMurtry JP, Thivierge MC, Davis TA (2004) Peripheral leptin administration alters hormone and metabolite levels in the young pig. *Comp Biochem Physiol* 138:17-25

Further references can be found on our website: www.mercodia.com

SUMMARY OF PROTOCOL SHEET**Merckodia Porcine Insulin ELISA**

Add Calibrators, controls* and samples	25 µL
Add enzyme conjugate 1X solution	100 µL
Incubate	2 hours at 18-25°C on a plate shaker (700-900 rpm)
Wash plate with wash buffer 1X solution	700 µL, 6 times
Add Substrate TMB	200 µL
Incubate	15 minutes at 18-25°C
Add Stop Solution	50 µL Shake for 5 seconds to ensure mixing
Measure A ₄₅₀	Evaluate results

*not provided

For full details see page 6